




SINGLE USE SURGICAL

PROFESSIONAL RANGE



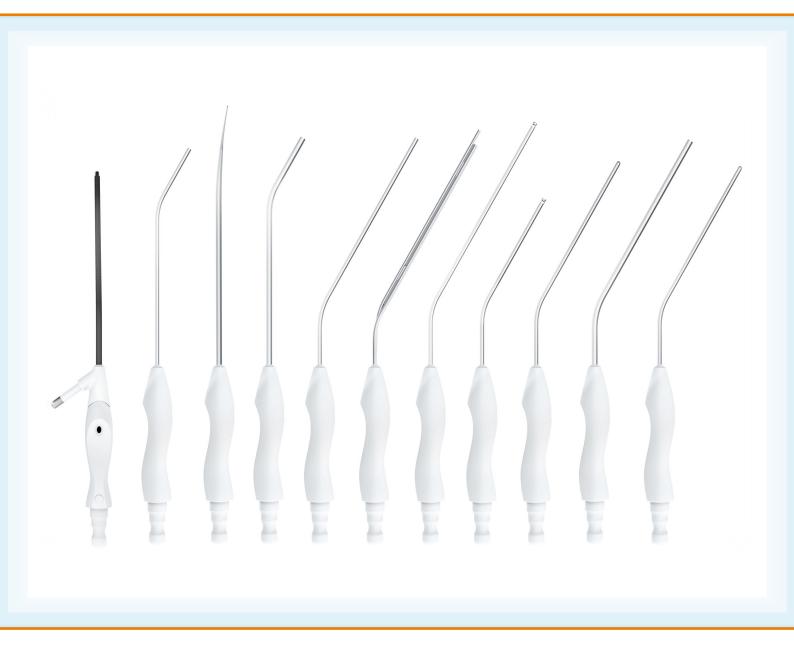

CONTENTS

SUCTIONS

CANNULAS

METAL INSTRUMENTS


OTOLOGY INSTRUMENTS

PULSE LAVAGE

| The Professional Range | 4 | OTOLOG | Y INSTRUMENTS |
|--------------------------------|-------|--------------------------------|---------------|
| Suctions | 5-27 | 5.5255 | PULSE LAVAGE |
| Suction Design | 5 | | 1 OLSE LAVAGE |
| Suction Size Chart | 6 | | |
| Polished Suctions | 7 | Sinus Trocar Cannula | 30 |
| Rounded Suctions | 8 | Spackman Cannula | 31 |
| Fenestrated Suctions | 9 | General Metal Instruments | 32-38 |
| Fenestrated & Rounded Suctions | 10 | Tischler-Morgan Biopsy Pun | ch 32 |
| Long Suctions | 11-12 | Cushings Artery Forceps | 33 |
| Teardrop Suctions | 13 | Jefferson Retractor | 33 |
| Fine Suctions | 14 | Mollinsons Retractor | 33 |
| Soft Tip Suction | 14 | Hudson Brace | 34 |
| Retractor Suction | 15 | Gigli Saw | 34 |
| Dissector Suction | 15 | Yasargil Retractor | 34 |
| Micro Suction | 16 | Fine Suture Forceps | 35 |
| Curved Micro Suction | 16 | Iris Scissors Straight & Curve | d 35 |
| Sinus Suctions | 17 | Kerrisons | 36 |
| Nasal Aspirating Dissector | 18 | Phlebectomy Instruments | 37 |
| Schuknechts | 18 | Mosquito Forceps | 37 |
| Curved Suction | 19 | Vein Hooks | 37 |
| Poole Tip Suction | 19 | Tracheal Dilator | 38 |
| Poole Suction | 19 | Tilley Nasal Forceps | 38 |
| Yankauer Suction | 20 | Thudichum Specula | 38 |
| Sheffield Suction | 20 | Otology Instruments | 39-42 |
| Laser Suctions | 21 | Myringotomy Knife | 39 |
| Rounded Fenestrated Suction | 22 | Myringotomy Kit | 39 |
| Weir Suction | 22 | Jobson Horne Probe | 40 |
| Insulated Suctions | 23 | Cerumen Hook | 40 |
| Suction Tubing | 23 | Rosen Inserter | 41 |
| Diathermy Suctions | 24-26 | Micro Ear Forceps | 41 |
| Diathermy Leads | 25 | Ear Specula | 42 |
| Sigmoidoscope Suction | 27 | Pulse Lavage | 43-44 |
| Aspiration Tube | 28 | Packaging | 45-46 |
| Aspiration Needle | 28 | Index by Product Reference | 47-48 |
| Cannula | 29-31 | Other Eakin Product Ranges | 51 |
| Tibbs Cannula | 29 | | |
| Cannula Straight | 29 | | |
| Cannula Curved | 29 | | 0.1.1 |
| Angled Cannula | 29 | | October 2021 |



THE PROFESSIONAL RANGE

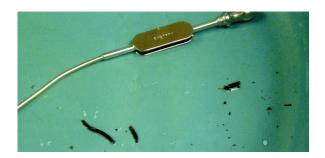

The problem: Debris found inside a reusable suction tube *after* cleaning and sterilization.



The solution: Single use suction tubes

GUARANTEED CLEAN

All cannulated instruments are manufactured in UK clean rooms, using processes that ensure internal lumens are guaranteed clean and sterile to eliminate the risk of contamination from manufacturing debris.

Every instrument is delivered to you in individual sterile pouches for immediate use in the sterile field.

CUSTOMER FOCUS

Meeting customer needs is at the heart of the company. Customer demand has driven the development of the extensive PROFESSIONAL range of suction tubes and other devices for many areas of surgery.

Eakin Surgical's friendly customer service team is always happy to work with you to find the right instrument to meet your requirements.

SUPERIOR QUALITY

Eakin Surgical strongly believes in value for money; it supplies the highest quality single use devices at competitive prices.

All products comply with industry standards as well as the exacting standards of professionals.

CONTACT US

Eakin Surgical Ltd Greypoint Cardiff Business Park Cardiff CF14 5WF United Kingdom

Telephone +44 (0)2920 767 800 Facsimile +44 (0)2920 767 801 Email info@eakinsurgical.com Website www.eakinsurgical.com

COST EFFECTIVE

When reusable instruments cannot be guaranteed clean, or are being lost, damaged and frequently replaced, single use instruments are a convenient cost effective option.





STANDARD SUCTION DESIGN

EXTENSIVE RANGE

The PROFESSIONAL range of suctions offers an extensive choice of length, lumen size and tip finish to enable you to find the right match for procedure, patient and surgeon.

Unless otherwise specified, all suctions have a 30° bend (a) to optimise visibility and this angle is fixed. The working length is measured from the bend to the tip (c).


ERGONOMIC HANDLE

A unique handle design provides optimum ergonomics with slim, smooth contours for comfort and improved visibility. Engineered to prevent blocking.

CONTROL HOLE

The suction control hole allows fine manipulation of suction with minimal noise. As standard the control hole is round but we also produce a range with a teardrop shaped control hole where more sensitive adjustment is needed (see page 13). There are also a small number of suctions that have no control hole in the handle, these will be indicated next to the listing.

| Ref. | Description |
|------|------------------------|
| а | Bend angle |
| b | Lumen/Diameter |
| С | Working length |
| d | Length before the bend |

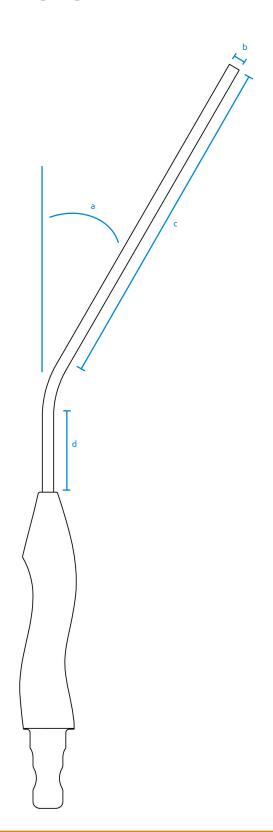

Suction Size & Tip Style Chart









| | Tip | | | | Polishe | ed | | | | | Rounded | | F | enestrate | ed | R Fe | ounded enestrate | & ed |
|---------------------|-----|---------|---------|---------|---------|---------|----------|------|-------|------|---------|------|------|-----------|------|---------|---------------------|---------|
| lumen | fg | 5 | 6 | 7 | 8 | 9 | 10 | 12 | 15 | 6 | 9 | 12 | 6 | 9 | 12 | 6 | 9 | 12 |
| ายก | mm | 1.7 | 2 | 2.3 | 2.7 | 3 | 3.3 | 4 | 5 | 2 | 3 | 4 | 2 | 3 | 4 | 2 | 3 | 4 |
| | 5/6 | 5fgN100 | 6fgN100 | 7fgN100 | 8fgN100 | 9fgN100 | 10fgN100 | | | | | S236 | S086 | S087 | S160 | | | |
| | 7 | | | | | | | | | S102 | S103 | | | | | | | |
| | 9 | 5fgN140 | 6fgN140 | 7fgN140 | 8fgN140 | 9fgN140 | 10fgN140 | S034 | S340 | S037 | S014 | S104 | S002 | S060 | | S091 | S092 | S093 |
| | 13 | | S077 | | | S076 | | S074 | S390 | S095 | S096 | S133 | S055 | S001 | S042 | S071 | S094 | S343 |
| | 15 | | | | | S005 | | | | | | | | | | | | |
| _ | 17 | | | | | | | S199 | | | | | | | | S120 | S075 | |
| Working length (cm) | 20 | | S031 | | | S024 | | S021 | S265* | S156 | S319 | S318 | | | | S168 | S169 | |
| ing le | 23 | | | | | | | S459 | | | | | | | | | | |
| ngth | 27 | | S020 | | | S011 | | | | S185 | S182 | | | | | | | |
| - (cm) | 34 | | | | | | | S460 | | | | | | | | | | |
| | 35 | | | | | S036 | | | | | S187 | S284 | | | | | | |
| | 40 | | | | | | | S109 | | | | | | | | | | |
| | 44 | | | | | | | S461 | | | | | | | | | | |
| | 45 | | S258 | | | S015 | | | | | S183 | | | | | | | |
| | 50 | | | | | | | S132 | | | | | | | | | | |
| | 55 | | | | | S111 | | S038 | | | | S184 | | | | | | |

Polished Tips

All suctions have smooth, polished tips.

For rounded and fenestrated tips, please see pages 8-10. See also our range of specialist suctions on pages 12-27.



| | Description | Length | Lumen | Ref. |
|---|-------------|--------|-------|----------|
| | Zoellner | 5.5cm | 5fg | 5fgN100 |
| DATE OF THE PARTY | | | 6fg | 6fgN100 |
| | | | 7fg | 7fgN100 |
| | Mastoid | 5.5cm | 8fg | 8fgN100 |
| | | | 9fg | 9fgN100 |
| | | | 10fg | 10fgN100 |
| | Frazier | 9cm | 5fg | 5fgN140 |
| | | | 6fg | 6fgN140 |
| | | | 7fg | 7fgN140 |
| | | | 8fg | 8fgN140 |
| | | | 9fg | 9fgN140 |
| | | | 10fg | 10fgN140 |
| | | | 12fg | S034 |
| | | | 15fg | S340 |
| | | 13cm | 6fg | S077 |
| | | | 9fg | S076 |
| | | | 12fg | S074 |
| | | | 15fg | S390 |
| | | 15cm | 9fg | S005 |
| | | 17cm | 12fg | S199 |
| | Magill | 3.5cm | 6fg | 6fgF140 |
| | | | 9fg | 9fgF140 |

Rounded Tips

A range of suction tubes with a smooth rounded tip for atraumatic tissue contact.

For rounded tip suctions with a working length of 20cm and over, please see page 11.



| Description | Length | Lumen | Ref. |
|-----------------|--------|-------|------|
| Rounded Suction | 5.5cm | 12fg | S236 |
| Rounded Frazier | 7cm | 6fg | S102 |
| | | 9fg | S103 |
| Rounded Frazier | 9cm | 6fg | S037 |
| | | 9fg | S014 |
| | | 12fg | S104 |
| Rounded Frazier | 13cm | 6fg | S095 |
| | | 9fg | S096 |
| | | 12fg | S133 |
| Rounded Magill | 3.5cm | 9fg | S046 |



Fenestrated Tips

A range of suctions with fenestrated tips to disperse aspiration and prevent suction trauma.



| Description | Length | Lumen | Ref. |
|---------------------|--------|-------|------|
| Fenestrated Suction | 5cm | 6fg | S086 |
| | | 9fg | S087 |
| | | 12fg | S160 |
| Fenestrated Suction | 9cm | 6fg | S002 |
| | | 9fg | S060 |
| Fenestrated Suction | 13cm | 6fg | S055 |
| | | 9fg | S001 |
| | | 12fg | S042 |

Fenestrated & Rounded Tips

Suctions that benefit from a tip that is both fenestrated and rounded.

For rounded, fenestrated suctions with a working length of 20cm and over, please see page 12.




| Description | Length | Lumen | Ref. |
|-----------------------------------|---|---|--|
| Fenestrated Rounded Frazier | 5cm | 6fg | S456 |
| | | 9fg | S457 |
| | | 12fg | S458 |
| | 9cm | 6fg | S091 |
| | | 9fg | S092 |
| | | 12fg | S093 |
| Fenestrated Rounded Frazier | 13cm | 6fg | S071 |
| i i dziei | | 9fg | S094 |
| | | 12fg | S343 |
| Fenestrated Rounded Suction | 17cm | 6fg | S120 |
| | | 9fg | S075 |
| | Fenestrated Rounded Frazier Fenestrated Rounded Frazier Fenestrated Rounded Frazier | Fenestrated Rounded Frazier 9cm Fenestrated Rounded Frazier Fenestrated Rounded Frazier 13cm Fenestrated Rounded Frazier | Fenestrated Rounded Frazier 5cm 6fg 9fg 12fg 6fg 9cm 9fg 12fg 12fg Fenestrated Rounded Frazier 13cm 6fg 9fg 12fg 12fg 12fg |



Long Suctions

A range of long suction tubes to optimize reach and visibility during ENT procedures.

- Brushed stainless steel prevents glare and reflection.
- Added benefit of rounded tips for extra atraumatic tissue contact.
- Strong high quality tube.

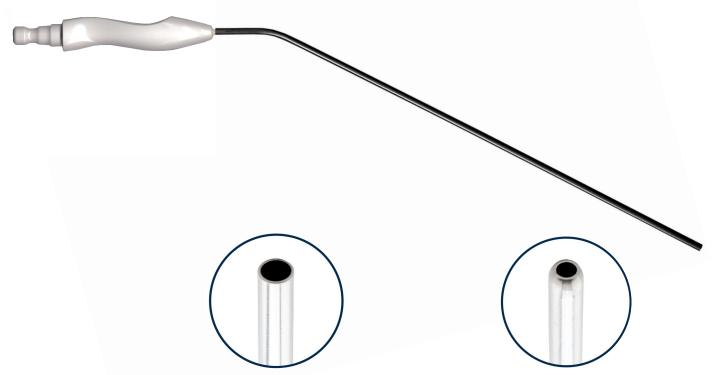

| Working | | Polishe | d Tip | | | Rounded Tip | |
|---------|-----------|-----------|------------|------------|----------|-------------|------------|
| Length | 6fg ~ 2mm | 9fg ~ 3mm | 12fg ~ 4mm | 15fg ~ 5mm | 6fg ~2mm | 9fg ~ 3mm | 12fg ~ 4mm |
| 20 cm | S031 | S024 | S021 | S265* | S156 | S319 | S318 |
| 23 cm | | | S459 | | | | |
| 27 cm | S020 | S011 | | | S185 | S182 | |
| 34 cm | | | S460 | | | | |
| 35 cm | | S036 | | | | S187 | S284 |
| 40 cm | | | S109 | | | | |
| 44 cm | | | S461 | | | | |
| 45 cm | S258 | S015 | | | | S183 | |
| 50 cm | | | S132 | | | | |
| 55 cm | | S111 | S038 | | | | S184 |

Made to order * Handle does not have a control hole

Long Suctions with 60° Bend

- ullet 60 $^\circ$ bend, provides a more acute angle than our standard suctions
- Brushed stainless steel prevents glare and reflection
- Strong high quality tube






| Working Length | Polished Tip | Rounded & F | enestrated Tip | |
|----------------|--------------|---------------|----------------|--|
| Working Length | 12fg ~ 4mm | 6fg ~ 2mm 9fg | | |
| 20 cm | | S168 | S169 | |
| 30 cm | S131 | | | |




Teardrop Control Suctions

- Teardrop shaped control hole for precise suction regulation
- Malleable tube gives the flexibility to adjust the angle by up to 50° for access to difficult to reach areas
- Available with a tapered tube to minimise blocking

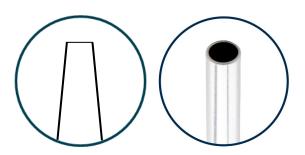



| Description | Length | Lumen | Ref. |
|--------------------------|--------|---------------------------|----------------------------------|
| Teardrop Micro Suction | 7cm | 18g | S439 |
| Teardrop Suction Tapered | 9cm | 6fg | S190 |
| Teardrop Suction Tapered | 13cm | 6fg | S193 |
| Teardrop Suction | 9cm | 7fg 8fg 9fg 10fg | \$452 \$453 \$189 \$454 |
| Teardrop Suction | 13cm | 12fg 9fg | S455 S192 |

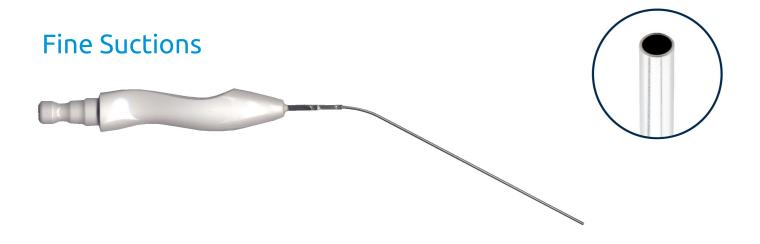

| Reference | Description | Length | Lumen |
|-----------|-----------------|--------|-------|
| S176 | Suction 3fg 8cm | 8cm | 3fg |
| S177 | Suction 5fg 8cm | 8cm | 5fg |



| Reference | Description | Length | Lumen |
|-----------|--------------------------|--------|-------|
| S321 | Soft Tip 6fg 9cm Suction | 9cm | 6fg |



Retractor Suction

A suction with a curved lip at its tip to allow simultaneous retraction and aspiration, available with teardrop control option.

The tube is flattened for easy access to the surgical site






| Reference | Description | Length |
|-----------|--|--------|
| S435 | Retractor Suction 12fg 12cm 60° | 12cm |
| S450 | Teardrop Retractor Suction 12fg 12cm 60° | 12cm |

Dissector Suction with Flattened Tip



| Reference | Description | Length |
|-----------|---|--------|
| S436 | Dissector Suction Flattened Tip 12fg | 13cm |
| S442 | Dissector Suction Flattened Tip 12fg 60 degrees | 10cm |

Made to order

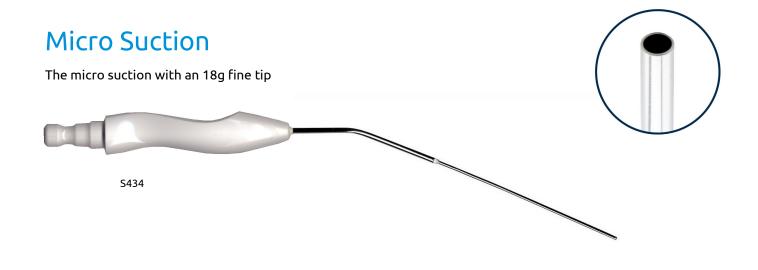

Curved Micro Suction

S446

Anatomically designed to access harder to reach areas. Perfect for procedures performed through the nose, the smooth rounded tip allows for gentle, atraumatic dissection.

30°



60°



| Reference | Description | Length | Lumen |
|-----------|--|--------|-------|
| S434 | Micro Suction 18g 7cm | 7cm | 18g |
| S446 | Curved Micro Suction 9fg 15cm 30° Bend Rounded Tip 60° Curve (Left) | 15cm | 9fg |
| S447 | Curved Micro Suction 9fg 15cm 30° Bend Rounded Tip 60° Curve (Right) | 15cm | 9fg |




Sinus Suctions

A range of suctions with curved tips for maximum visibility and access in difficult to reach areas





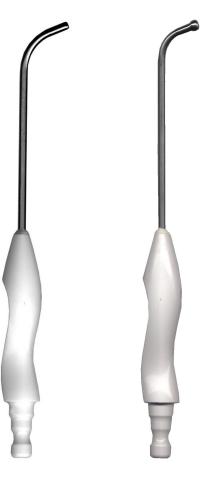



| | Curve Rounded Tip | | Curve Ball End | | Curve Rounded Tip | | Curve Ball End | |
|----------------------|-------------------|------------|----------------|------------|-------------------|------------|----------------|------------|
| | 70 |)° | 70° | 65° | 9 | 0° | 9 | 0° |
| | 9fg ~ 3mm | 12fg ~ 4mm | 9fg ~ 3mm | 12fg ~ 4mm | 9fg ~3mm | 12fg ~ 4mm | 9fg ~ 3mm | 12fg ~ 4mm |
| Upward Curve | | | S467 | S468 | | | | |
| With Control Hole | S012 | S013 | S228 | S229 | S040 | S041 | S230 | S231 |
| No Control Hole | | S163 | | | | | | |


Nasal Aspirating Dissector

The flattened curved tip of this instrument makes it ideal for access and elevation in the nasal cavity. A small fenestration on both sides of the tip allows simultaneous suction






| Reference | Description | Length |
|-----------|----------------------------|--------|
| S146 | Nasal Aspirating Dissector | 12cm |

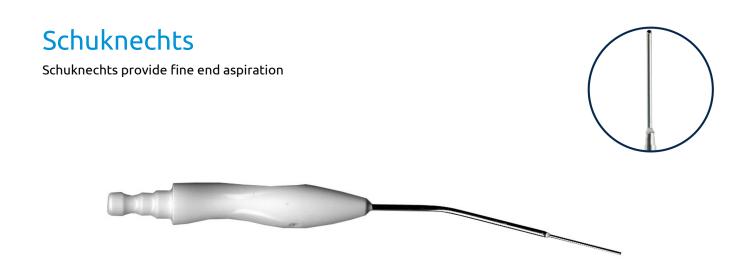

| Reference | Description | Length | Lumen |
|-----------|---------------------|--------|-------|
| S007 | Schuknecht 20g | 6cm | 20g |
| S008 | Schuknecht 24g | 6cm | 24g |
| S072 | Schuknecht 18g | 6cm | 18g |
| S157 | Long Schuknecht 20g | 13cm | 20g |



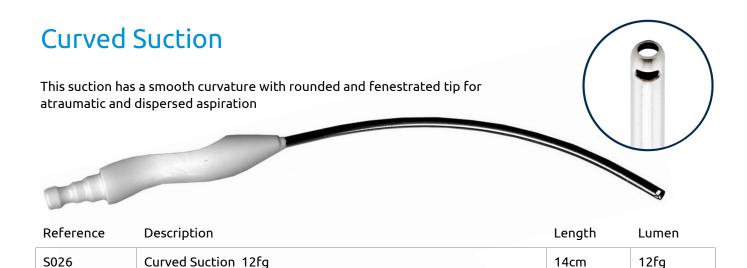




| Reference | Description | Length | Lumen |
|-----------|---|--------|-------|
| S440 | Poole Tip Suction 12fg 13cm Rounded 4x4x2mm Fenestrations | 13cm | 12fg |

Poole Suction

The Poole suction has multiple holes around the sheath for dispersed suction and reduced risk of blockage. It also offers the option for finer suction with a simple bayonet twist action



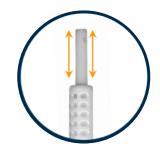

| Reference | Description | Length | Lumen |
|-----------|---------------|--------------------|-------|
| S010 | Poole Suction | 26cm (when closed) | 5mm |

Yankauer Suction



Available in ball and fenestrated tip style variants and ensures reliability throughout the entire procedure



| Reference | Description |
|-----------|----------------------------------|
| GC1233K3W | Yankauer Suction Ball Tip |
| GC1215A3W | Yankauer Suction Fenestrated Tip |

 ${\it Yankauer Suctions are under distribution agreement of GCMEDICA\ Enterprise\ Ltd.}$

Sheffield Suction

Suction with novel clip to allow attachment to a laryngoscope during endolaryngeal laser surgery.



Enables effective smoke plume evacuation and is detachable in the instance of a bleed



| Reference | Description | Length | Lumen |
|-----------|---|--------|-------|
| S240 | Sheffield Suction 9fg 15cm 90° Bend with Clip | 15cm | 9fg |



Laser Suctions

Laser suctions combine two instruments in one; a 0.9mm fibre optic laser port bonded to the suction tube allows simultaneous aspiration with cutting and coagulation

For use with 600 micron core laser fibre with 850 micron sheath



Extended Laser Channel



Flush Finish Laser Channel



Extended Angulated Laser Channel

| Description | Length | Lumen | Ref. |
|--|--------|-------|------|
| Laser Suction 6fg 9cm with Extended Laser Channel | 9cm | 6fg | S110 |
| Laser Suction with Flush Finish | 9cm | 6fg | S165 |
| | 20cm | 6fg | S139 |
| Fairley Laser Suction 12cm with Extended Angulated Tip | 12cm | 6fg | S285 |

Rounded Suction with 4 Fenestrations

This ultra delicate suction has 4 fenestrations that are 1mm in diameter at the tip





| Reference | Description | Length | Lumen |
|-----------|---|--------|-------|
| S147 | Rounded Suction 12fg 9cm 4x1mm Fenestrations | 9cm | 12fg |
| S445 | Rounded Suction 12fg 30cm 4x1mm Fenestrations 60° | 30cm | 12fg |

Weir Suction

A suction with a flattened, semi-sharp open end to combine minor dissection with aspiration






| Reference | Description | Length |
|-----------|-------------------|--------|
| S009 | Weir Suction 12fg | 13cm |



Insulated Suctions

A range of suctions with protective insulation to allow use alongside diathermy



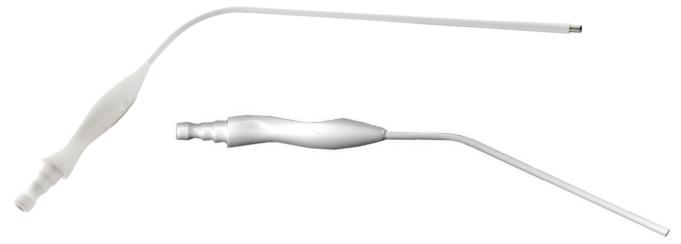

| Reference | Description | Length | Lumen |
|-----------|--------------------------------|--------|-------|
| S432 | Insulated Suction 9fg 17cm 90° | 17cm | 9fg |
| S433 | Insulated Suction 6fg 17cm 90° | 17cm | 6fg |
| S025 | Insulated Frazier 9fg 9cm | 9cm | 9fg |

Made to order

Flexible Connection Tubing

Made from silicone for flexibility and resistance to kinking



| Reference | Description | Includes Connector |
|-----------|--------------------------------------|--------------------|
| S066 | Flexible Connection Tubing 8mm 50cm | No |
| S099 | Flexible Connection Tubing 8mm 300cm | Yes |
| S101 | Flexible Connection Tubing 8mm 150cm | Yes |

Diathermy Suctions

As a 2 in 1 instrument, diathermy suctions create a clearer view and a less invasive procedure by minimising the number of instruments in the field of vision.

The suction helps remove harmful smoke from theatre. The 4mm pin connects to diathermy lead (see opposite page for reusable leads).





| Working | Polished Tip | | Rounded Tip | Fenestrated Tip | Rounded & Fenestrated Tip |
|---------|--------------|-----------|-------------|-----------------|------------------------------|
| Length | 9fg ~3mm | 12fg ~4mm | 9fg ~ 3mm | 9fg ~3mm | 9fg ~ 3mm |
| 9 cm | S117 | | S136 | | |
| 13 cm | S116 | | | | S119 |
| 19 cm | S195 | | | | |
| 20 cm | S115 | | | | |
| 22 cm | | | | S113 | |
| 27 cm | S114 | | | | |





Featuring a 70° bend and a standard polished tip.



Features on-tip smoke evacuation as recommended by HSE

Diathermy Leads

Competitively priced reusable leads for use with Single Use Surgical 'Professional' Diathermy suctions



| Reference | Description | Length | Lumen |
|-----------|-----------------------------------|--------|-------|
| S438 | Diathermy Suction 9fg 17cm 70° | 17cm | 9fg |
| 1011-35 | Diathermy lead 4mm pin 4mm socket | 3m | n/a |
| 1012-35 | Diathermy lead 8mm pin 4mm socket | 3m | n/a |

Diathermy Suctions

Anatomically designed to access harder to reach areas often associated with Submucosal Diathermy (SMD).

Combined diathermy and suction for targeted and effective smoke extraction

- Built in smoke extraction

- Ø3mm—6fg 10cm
- 40° bend
- 4mm connection pin
- Latex & DEHP free



I have used the Diathermy Abbey Needle with Suction, and I like it. It is quicker and less traumatic to use. It feels more robust than the usual needle, so you can push the turbinate laterally as well, creating a bit more airway and avoiding accidental diathermy of the septum. This reduces the chance of complications such as intra-nasal adhesion occurring. Simultaneous clearance of the smoke means you don't need to alternate between the diathermy needle and the Zoellner sucker, which can often cause minor bleeding. It allows good vision throughout the procedure.

David Hill, Consultant ENT Surgeon

| Reference | Description | Length | Lumen |
|-----------|-------------------------------------|--------|-------|
| S448 | Diathermy Abbey Needle with Suction | 10cm | 6fg |



Sigmoidoscope Suction

A rigid sigmoidoscope suction without a suction control hole





| Reference | Description | Length | Lumen |
|-----------|---------------------------------|--------|-------|
| S043 | Sigmoidoscope Suction 12fg 44cm | 44cm | 12fg |

Aspiration Needles

Aspiration needles with a sharp bevelled tip for safe and controlled penetration when aspirating fluids from cysts and organs. With Luer Lock. Please note, product colour may vary from that pictured.







| Reference | Description | Length | |
|-----------|----------------------------|--------|--|
| S019 | Aspiration Needle 3mm 32cm | 32cm | |
| S080 | Aspiration Needle 5mm 32cm | 32cm | |

Aspiration Tube

Atraumatic blunt ended tube with Luer Lock





| Reference | Description | Length | Lumen |
|-----------|--------------------------------|--------|-------|
| S112 | Blunt Aspiration Tube 4mm 33cm | 33cm | 4mm |



Tibbs Cannula

The Tibbs hybrid cannula combines traditional olive and cone tips




| Reference | Description | Length |
|-----------|------------------------|--------|
| S044 | Tibbs Cannula Set of 3 | 7cm |



| Reference | Description | Length | Lumen |
|-----------|-----------------------------|--------|-------|
| S214 | 6fg Cannula Straight | 7cm | 6fg |
| S213 | 6fg Cannula Curved Set of 2 | 7cm | 6fg |
| S028 | Angled Cannula 19g | 2.5cm | 19g |

Sinus Trocar Cannula

The Sinus Trocar and Cannula - also known as an antrum cannula - consists of a removable cannula over a 11.2cm trocar to allow washout of the sinuses.

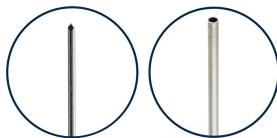



| Reference | Description |
|---------------|-------------|
| I CI CI CIICC | Description |

| S081 Sinus Trocar Cannula | |
|---------------------------|--|
|---------------------------|--|



Spackman Cannula

The single use Spackman Cannula combats concerns over cleanliness with the added benefit of an easy to use self-locking clamp. Available with a fixed or repositionable acorn.

For hydrotubation dye tests and uterine manipulation during laparoscopy.





Fixed Acorn

Moveable Acorn




| Reference | Description | Length | Lumen |
|-----------|-------------------------------|--------|-------|
| S029 | Spackman Cannula | 32cm | 12fg |
| S030 | Spackman Cannula Moving Acorn | 32cm | 12fg |

Tischler-Morgan Biopsy Punches

The biopsy punches come with either a straight tip (pictured right) or a tip that is angled down (pictured below) to allow surgical staff to select the best option for their specific patient and procedure.

- Single use eliminates patient cross infection risk, biopsy sample cross contamination and guarantees quality
- Tip locking system safely secures sample
- Ergonomic hand-polished edges
- Quality & functionality matches reusables
- Bite: 7mm x 3.5mm
- Working Length: 200mm



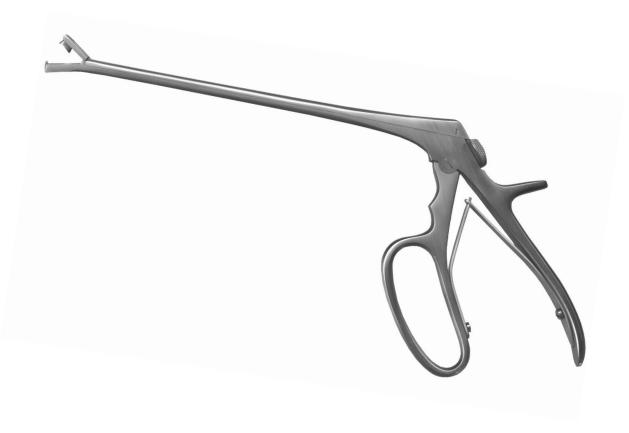

| Reference | Description | |
|-----------|--|--|
| S263 | Tischler-Morgan Biopsy Punch 200mm, Straight Tip | |
| S320 | Tischler-Morgan Biopsy Punch 200mm, Tip Down | |





Reference Description

S205 Cushings Artery Forceps Curved, Half Serrated Jaws, 150mm

Made to order



Reference Description

S206 Jefferson Retractor, Self-Retaining, Curved, Blunt Pointed Blades, 125mm




Reference Description

S207 Mollinsons Retractor, Self-Retaining, Curved, 4/4 Teeth, Adult, 165mm

Hudson Brace, Gigli Saw and Yasargil Retractor



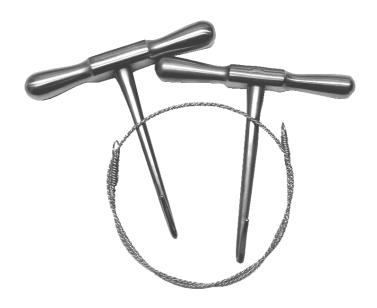

| Reference | Description | |
|-----------|---|--|
| S211 | Hudson Brace, Standard Size with Universal Fitting, 265mm | |
| S212 | Gigli Saw 300mm and Handles, Pair, 55mm with Guide | |
| S419 | Yasargil Retractor (not pictured) | |



Fine Suture Forceps



| Reference | Description |
|-----------|---------------------|
| S257 | Fine Suture Forceps |

Made to order


Iris Scissors Straight & Curved



| Reference | Description |
|-----------|------------------------|
| S255 | Iris Scissors Straight |
| S256 | Iris Scissors Curved |

Kerrisons

- Designed for single patient use to eliminate cross infection and guarantee a sharp cut
- Ergonomic hand-polished edges
- Quality & functionality matches reusables

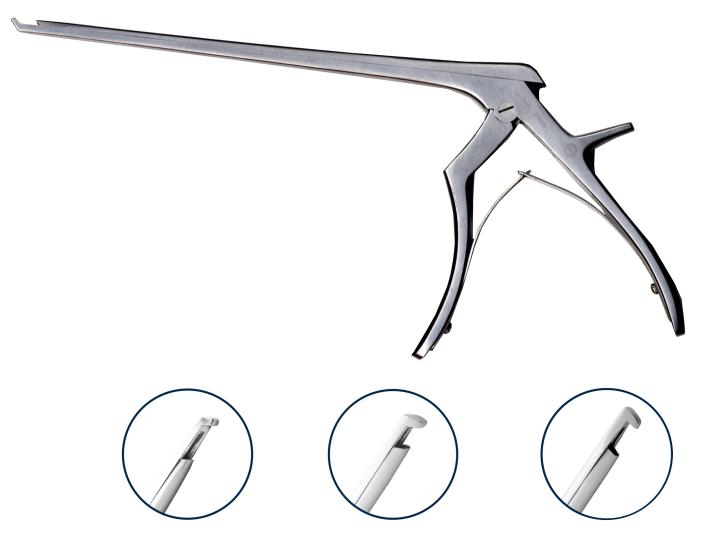

| Bite | 45° Degree Upward Cutting (Delicate Use) | 45° Degree Upward Cutting | 90° Degree |
|------|---|---------------------------|------------|
| 1mm | S161 | S125 | |
| 2mm | S261 | S126 | S367 |
| 3mm | S262 | S127 | |
| 4mm | | S137 | |
| 5mm | | S138 | |



Vein Hooks

A cost effective alternative to reusable equivalents. Made from high grade stainless steel they have a strong Oesch style hook which can withstand rigorous handling, yet is fine enough for insertion through small incisions





We have recently started using phlebectomy hooks and mosquito clips from Eakin Surgical Ltd at all of our clinics. The first thing we noticed was the very good quality. For disposable surgical instruments, the feel and performance of the hooks and the clips is excellent. As the quality of our work depends upon the quality of the instruments that we use, we are delighted to have found Eakin Surgical Ltd.

Professor Mark Whiteley, Venous Surgeon and owner of the Whiteley Clinic

| Reference | Description |
|-----------|-------------------------------------|
| S090 | Medium Right Hand Vein Hook (No. 2) |
| S106 | Small Right Hand Vein Hook (No. 3) |
| S134 | Medium Left Hand Vein Hook (No. 2) |
| S135 | Small Left Hand Vein Hook (No. 3) |

Mosquito Forceps

Single use mosquito forceps provide convenience and consistent quality



Made to order

| Reference | Description |
|-----------|----------------------------------|
| S098 | Mosquito Forceps Curved Set of 2 |



| Reference | Description |
|-----------|-------------|
| | |

S181 Tracheal Dilator

Made to order





| Reference | Description |
|-----------|-------------|
| Reference | Description |

| S180 | Tilley Nasal Forceps Short |
|------|----------------------------|
| S380 | Tilley Nasal Forceps Long |

Thudichum Specula



| Reference | Description |
|-----------|--------------------------|
| S179 | Thudichum Specula Size 5 |
| S238 | Thudichum Specula Size 4 |
| S239 | Thudichum Specula Size 3 |



Myringotomy Kit

These all in one procedure packs contain sterilised, easy to use disposable instruments that are packed in a sterile peel pouch to save on cleaning and assembly costs



| Reference | Description |
|-----------|---|
| S301 | Myringotomy Kit 1 (Micro Forcep, Knife, Jobson Horne Probe, Rosen Inserter, Ear Specula 3/4/5/6mm, Zoellner Suction Tube, 18g Zoellner Tip, 22g Zoellner Tip) |

Myringotomy Knife

Lancet used to pierce the tympanic membrane. Bendable shaft for ease, versatility and to prevent visual obstruction. Handle is shaped to aid orientation of the blade.






| Reference | Description | Total length | Handle length |
|-----------|-------------------|--------------|---------------|
| S310 | Myringotomy Knife | 157mm | 95mm |

Jobson Horne Probe

Designed to allow gentle and easy extraction of wax or foreign bodies from the ear canal. The probe has a Curette loop at one end and a threaded section at the other that can be used to hold wrapped cotton wool.

The handle section in the centre features flattened sides to aid grip and prevent unintentional twisting. Made from carbon fibre filled nylon for strength and flexibility.

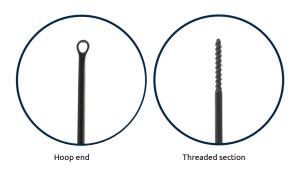




| Reference | Description | Length | Loop |
|-----------|----------------------------------|--------|------|
| S303 | Jobson Horne Probe | 150mm | 3mm |
| S311 | Jobson Horne Probe (Non-Sterile) | 150mm | 3mm |

Cerumen Hook

Disposable cerumen hook for extracting hard wax and other obstacles from the ear canal. The hook shape has been designed to offer good leverage whilst the flexible carbon fibre helps reduce trauma to the ear.



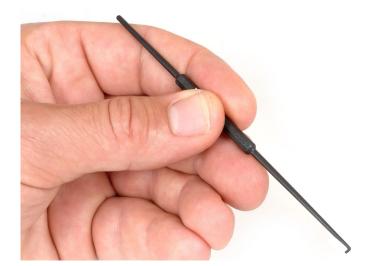

| Reference | Description | Length | Hook Length |
|-----------|----------------------------|--------|-------------|
| S307 | Cerumen Hook | 135mm | 2mm |
| S312 | Cerumen Hook (Non-Sterile) | 135mm | 2mm |



Rosen Inserter

Developed with surgeons to aid grommet insertion. The geometry also proves favourable for removal of objects in the ear.

Manufactured from carbon fibre filled nylon for strength and flexibility.



| Reference | Description | Length |
|-----------|------------------------------|--------|
| S308 | Rosen Inserter | 143mm |
| S317 | Rosen Inserter (Non-Sterile) | 143mm |

Made to order

Micro Ear Forceps

High quality micro ear forceps—also known as crocodile forceps—that effectively remove foreign bodies or are used to insert grommets. Feature serrated jaws to easily grip particles and fine shaft for access to the inner ear.



| Reference | Description | Shaft length |
|-----------|----------------------------|--------------|
| S221 | Micro Ear Forceps 3mm jaws | 85mm |

Ear Specula

Advanced design provides maximum patient comfort and superior operational vision

Matt finish and internal ridges reduce glare and reflection from the microscope light onto the internal surface of the specula. Whilst the end is slightly oval to conform to the ear canal

Soft and flexible material to aid patient comfort

Serrated lateral end aids finger tip grip



| Reference | Description |
|-----------|-------------------------------|
| S304 | Ear Specula 4mm |
| S305 | Ear Specula 5mm |
| S306 | Ear Specula 6mm |
| S309 | Ear Specula 3mm |
| S313 | Ear Specula 3mm (Non-Sterile) |
| S314 | Ear Specula 4mm (Non-Sterile) |
| S315 | Ear Specula 5mm (Non-Sterile) |
| S316 | Ear Specula 6mm (Non-Sterile) |

Tumarkin Ear Specula

Metal Tumarkin Ear Specula.

Serrated lateral end aids finger tip grip.




| Reference | Description |
|-----------|--------------------------------------|
| S322 | Tumarkin Ear Specula Size 0.5, 4.5mm |



Pulse Lavage Kit

A disposable pulse lavage system with a built-in battery, for controlled pulsatile irrigation and suction


- Built-in battery pack easily removed for safe disposal**; cutting the weight and cost of clinical waste
- Long coaxial tip provides extra reach
- 30mm cone tip
- Ergonomic handle

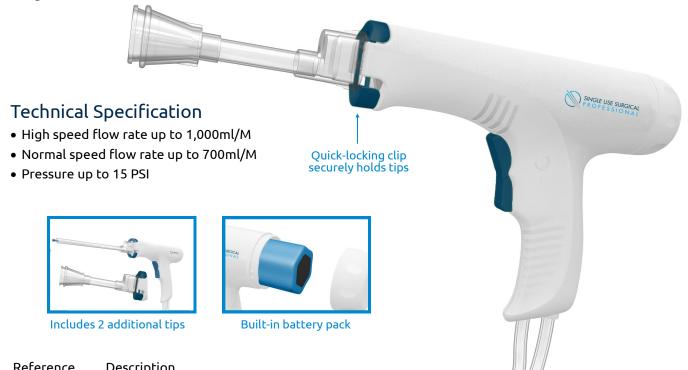

| Reference | Description | | |
|-----------------|---|--|--|
| S353 | Pulse Lavage Kit 1 includes Pulse Lavage, 30mm Cone Tip, Long Coaxial Canal (Femoral) Tip | | |
| Accessories - s | Accessories - see next page for images | | |
| BT 04 Z | Coaxial Canal (Femoral) Brush Tip | | |
| CT 02 Z | Coaxial Canal (Femoral) Tip | | |
| CTL 02 Z | Long Coaxial Canal (Femoral) Tip | | |
| LT 01 Z | 75mm Cone Tip | | |
| MT 01 Z | 50mm Cone Tip | | |
| ST 01 Z | 30mm Cone Tip | | |
| RT23012Z | Hip Revision Tip | | |

^{**} Ensure batteries are disposed of safely and in accordance with all current national and international waste disposal directives

Pulse Lavage - Accessories

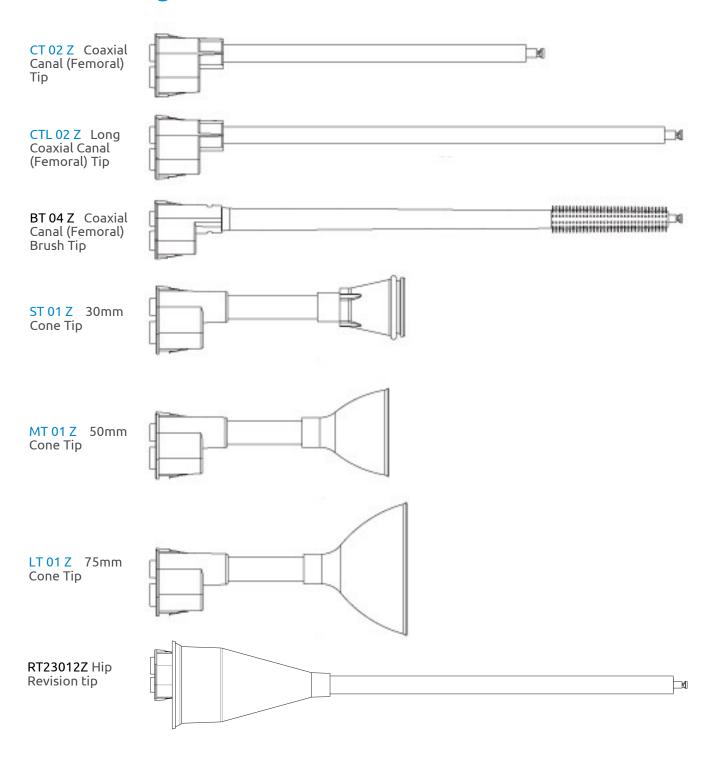

Pulse Lavage Tip accessories are under distribution agreement of Dongguan Kaiser Technology Co Ltd.



Packaging

STERILE POUCHES

All sterile instruments are supplied in individual peelable pouches to protect them during storage until needed.

DISPENSING BOXES

Each PROFESSIONAL box has a removable lid for convenient dispensing of sterile pouches even when boxes are stacked. The front label on the boxes has a barcode for quick identification and stock management.

CARE FOR THE ENVIRONMENT

Eakin Surgical strive to protect the environment through responsible and sustainable operations. All cardboard packaging is recyclable. For information on recycling of metal instruments please ask your representative

Packaging Symbol Key

| Manufacturer | STERILEEO | Sterilized using ethylene oxide | NON | Non– sterile |
|---|------------|---|--------------|--|
| REF Catalogue number | STERNIZE | Do not resterilize | * | Keep dry |
| LOT Batch code | | Single barrier system with protective packaging outside | UDI | Unique Device Identifier |
| Do not use if package is damaged and consult instructions for use | []i | Consult instructions for use or consult electronic instructions for use | | Latex Free |
| Do not re-use | MD | Medical Device | -10°C | Temperature limit |
| Caution | \square | Use-by-date | 촟 | Keep away from sunlight |
| Date of Manufacture | EC REP | Authorized representative in the European Community/European Union | | Refer to instruction manual/ booklet |
| QTY Quantity of content | \bigcirc | Single sterile barrier system | ※ | Protect from heat and radioactive sources |
| Indicates a recovery/ recyclable package or pack- age material | ALKALINE | Alkaline Battery (recyclable) | Ā | Recycle in accordance with The Waste Electrical and Electronic Equipment |
| R CAUTION: US Federal law restricts this device to sale by or on the order of the physician | | | | Directive (WEEE Directive) European Community directive 2012/19/EU |



CE mark including SGS (notified body) identification number. Product conforms with European Requirements



INDEX BY PRODUCT REFERENCE

| Reference | Description | Page | Reference | Description | Page |
|--------------|--|----------|-----------|--|------|
| 1011-35 | 3m Diathermy lead 4mm pin, 4mm socket. | 25 | S055 | Fenestrated Suction 6fg 13cm | 9 |
| 1012-35 | 3m Diathermy lead, 8mm pin, 4mm socket | 25 | S060 | Fenestrated Frazier 9fg 9cm | 9 |
| 10fgN100 | Mastoid Suction 10fg | 7 | S066 | Flexible Connection Tubing 8mm 50cm | 23 |
| 10fgN140 | Frazier Suction 10fg 9cm | 7 | S071 | Fenestrated Rounded Frazier 6fg 13cm | 10 |
| 5fgN100 | Zoellner Suction 5fg | 7 | 5072 | Schuknecht 18g 6cm | 18 |
| 5fgN140 | Fine Frazier Suction 5fg 9cm | 7 | 5074 | Suction 12fg 13cm | 7 |
| 6fgF140 | Fine Magill 6fg | 7 | S075 | Fenestrated Rounded Suction 9fg 17cm | |
| 6fgN100 | Zoellner 6fg | 7 | | - | 10 |
| 6fgN140 | Fine Frazier 6fg 9cm | 7 | S076 | Suction 9fg 13cm | 7 |
| 7fgN100 | Zoellner Suction 7fg | 7 | S077 | Suction 6fg 13cm | 7 |
| 7fgN140 | Fine Frazier Suction 7fg 9cm | 7 | S080 | Aspiration Needle 5mm 32cm | 28 |
| 8fgN100 | Mastoid Suction 8fg | 7 | S081 | Sinus Trocar/Cannula | 30 |
| 8fgN140 | Frazier Suction 8fg 9cm | 7 | S086 | Fenestrated Suction 6fg 5cm | 9 |
| 9fgF140 | Magill Suction 9fg | 7 | S087 | Fenestrated Suction 9fg 5cm | 9 |
| 9fgN100 | Mastoid Suction 9fg | 7 | S090 | Medium Right Hand Vein Hook (No. 2) | 37 |
| 9fgN140 | Frazier Suction 9fg 9cm | 7 | S091 | Fenestrated Rounded Frazier 6fg 9cm | 10 |
| BT 04 Z | Coaxial Canal (Femoral) Brush Tip | 44 | S092 | Fenestrated Rounded Frazier 9fg 9cm | 10 |
| CT 02 Z | Coaxial Canal (Femoral) Tip | 44 | S093 | Fenestrated Rounded Frazier 12fg 9cm | 10 |
| CTL 02 Z | Long Coaxial Canal (Femoral) Tip | 44 | S094 | Fenestrated Rounded Frazier 9fg 13cm | 10 |
| GC1215A3W | Yankauer Suction Ball Tip | 20 | S095 | Rounded Frazier 6fg 13cm | 8 |
| GC1233K3W | Yankauer Suction Fenestrated Tip | 20 | S096 | Rounded Frazier 9fg 13cm | 8 |
| LT 01 Z | 75mm Cone Tip | 44 | S098 | Mosquito Forceps Curved Set of 2 | 37 |
| MT 01 Z | 50mm Cone Tip | 44 | S099 | Flexible Connection Tubing 8mm 300cm | 23 |
| RT23012Z | Hip Revision Tip | 44 | S101 | Flexible Connection Tubing 8mm 150cm | 23 |
| ST 01 Z | 30mm Cone Tip | 44 | S102 | Rounded Frazier 6fg 7cm | 8 |
| S001 | Fenestrated Suction 9fg 13cm | 9 | S103 | Rounded Frazier 9fg 7cm | 8 |
| 5002 | Fenestrated Suction 6fg 9cm | 9 | S104 | Rounded Frazier 12fg 9cm | 8 |
| S005 | Suction 9fg 15cm | 7 | S106 | Small Right Hand Vein Hook (No. 3) | 37 |
| S007 | Schuknecht 20g 6cm | 18 | S100 | | |
| S008 | Schuknecht 24g 6cm | 18 | | Long Suction 12fg 40cm | 11 |
| S009 | Weir Suction 12fg | 22 | S110 | Laser Suction 6fg 9cm | 21 |
| S010 | Poole Suction | 19 | S111 | Long Suction 9fg 55cm | 11 |
| S011 | Long Suction 9fg 27cm | 11 | S112 | Blunt Aspiration Tube 4mm 33cm | 28 |
| S012 | Sinus 9fg | 17 | S113 | Diathermy Suction 9fg 22cm Fenestrated | 24 |
| S012 S013 | Sinus 12fq | 17 | S114 | Diathermy Suction 9fg 27cm | 24 |
| S014 | Rounded Frazier 9fg 9cm | 8 | S115 | Diathermy Suction 9fg 20cm | 24 |
| S014 S015 | Long Suction 9fg 45cm | 11 | S116 | Diathermy Suction 9fg 13cm | 24 |
| S019 | Aspiration Needle 3mm 32cm | 28 | S117 | Diathermy Suction 9fg 9cm Frazier | 24 |
| S020 | Long Suction 6fg 27cm | 11 | S119 | Diathermy Suction 9fg 13cm Fenestrated Rounded | 24 |
| | | | S120 | Fenestrated Rounded Suction 6fg 17cm | 10 |
| S021 S024 | Long Suction 12fg 20cm Long Suction 9fg 20cm | 11 | S125 | Kerrison Punch 1mm x 200mm @45° | 36 |
| S024 S025 | 3 3 | 11 23 | S126 | Kerrison Punch 2mm x 200mm @45° | 36 |
| | Insulated Frazier 9fg 9cm | | S127 | Kerrison Punch 3mm x 200mm @45° | 36 |
| S026 | Curved Suction 12fg | 19 | S131 | Long Suction 12fg 30cm @60° | 12 |
| 5028 | Angled Cannula 19g | 29 | S132 | Long Suction 12fg 50cm | 11 |
| 5029 | Spackman Cannula Maying Asses | 31 | S133 | Rounded Suction 12fg 13cm | 8 |
| 5030 | Spackman Cannula Moving Acorn | 31 | S134 | Medium Left Hand Vein Hook (No. 2) | 37 |
| 5031 | Long Suction 6fg 20cm | 11 | S135 | Small Left Hand Vein Hook (No. 3) | 37 |
| 5034 | Large Frazier 12fg 9cm | 7 | S136 | Diathermy Suction 9fg 9cm Rounded | 24 |
| 5036 | Long Suction 9fg 35cm | 11 | S137 | Kerrison Punch 4mm x 200mm @45° | 36 |
| 5037 | Rounded Fine Frazier 6fg 9cm | 8 | | | |
| 5038 | Long Suction 12fg 55cm | 11 | S138 | Kerrison Punch 5mm x 200mm @45° | 36 |
| 5040 | Sinus 9fg Small Bend | 17 | S139 | Long Laser Suction 6fg 20cm Flush Finish | 21 |
| 5041 | Sinus 12fg Small Bend | 17 | S146 | Nasal Aspirating Dissector | 18 |
| 5042 | Fenestrated Suction 12fg 13cm | 9 | S147 | Rounded Suction 12fg 9cm 4x1mm Fenestrations | 22 |
| S043 | Sigmoidoscope Suction 12fg 44cm | 27 | S156 | Rounded Long Suction 6fg 20cm | 11 |
| S044 | Tibbs Cannula Set of 3 | 29 | S157 | Long Schuknecht 20g 12.5cm | 18 |
| S046 | Rounded Magill 9fg | 8 | S160 | Fenestrated Suction 12fg 5cm | 9 |

INDEX BY PRODUCT REFERENCE

| Reference | Description | Page |
|-----------|--|------|
| S161 | Kerrison Punch Fine 1mm x 200mm for Delicate Use | 36 |
| S163 | Sinus 12fg No Control Hole | 17 |
| S165 | Laser Suction 6fg 9cm Flush Finish | 21 |
| S168 | Fenestrated Rounded Suction 6fg 20cm @60° | 12 |
| S169 | Fenestrated Rounded Suction 9fg 20cm @60° | 12 |
| S176 | Suction 3fg 8cm | 14 |
| S177 | Suction 5fg 8cm | 14 |
| S179 | Thudichum Specula Size 5 | 38 |
| S180 | Tilley Nasal Forceps Short | 38 |
| S181 | Tracheal Dilator | 38 |
| S182 | Rounded Long Suction 9fg 27cm | 11 |
| S183 | Rounded Long Suction 9fg 45cm | 11 |
| S184 | Rounded Long Suction 12fg 55cm | 11 |
| S185 | Rounded Long Suction 6fg 27cm | 11 |
| S187 | Rounded Long Suction 9fg 35cm | 11 |
| S189 | Teardrop Suction 9fg 9cm | 13 |
| S190 | Teardrop Suction 6fg 9cm Tapered | 13 |
| S192 | Teardrop Suction 9fg 13cm | 13 |
| S193 | Teardrop Suction 6fg 13cm Tapered | 13 |
| S195 | Diathermy Suction 12fg 19cm | 24 |
| S199 | Suction 12fg 17cm | 7 |
| S205 | Cushings Artery Forceps Curved, Half Serrated Jaws, 150mm | 33 |
| S206 | Jefferson Retractor, Self-Retaining, Curved, Blunt Pointed Blades, 125mm | 33 |
| S207 | Mollinsons Retractor, Self-Retaining, Curved, 4/4 Teeth, Adult, 165mm | 33 |
| S211 | Hudson Brace, Standard Size with Universal Fitting, 265mm | 34 |
| S212 | Gigli Saw 300mm and Handles, Pair, 55mm with Guide | 34 |
| S213 | 6fg Cannula Curved Set of 2 | 29 |
| S214 | 6fg Cannula Straight | 29 |
| S221 | Micro Ear Forceps 3mm | 41 |
| S228 | Ball End Sinus 9fg | 17 |
| S229 | Ball End Sinus 12fg | 17 |
| S230 | Ball End Sinus Small Bend 9fg | 17 |
| S231 | Ball End Sinus Small Bend 12fg | 17 |
| S236 | Rounded Suction 12fg 5.5cm | 8 |
| S238 | Thudichum Specula Size 4 | 38 |
| S239 | Thudichum Specula Size 3 | 38 |
| S240 | Sheffield Suction 9fg 15cm 90° Bend with Clip | 20 |
| S255 | Iris Scissors Straight | 35 |
| S256 | Iris Scissors Curved | 35 |
| S257 | Fine Suture Forceps | 35 |
| S258 | Long Suction 6fg 45cm | 11 |
| S261 | Kerrison Punch Fine 2mm x 200mm for Delicate Use | 36 |
| S262 | Kerrison Punch Fine 3mm x 200mm for Delicate Use | 36 |
| S263 | Tischler-Morgan Biopsy Punch 200mm, Straight Tip | 32 |
| S265 | Long Suction 15fg 20cm No Control Hole | 11 |
| S284 | Rounded Long Suction 12fg 35cm | 11 |

| Reference | Description | Page |
|-----------|--|------|
| S285 | Fairley Laser Suction 12cm with Extended Angulated Tip | 21 |
| S301 | Myringotomy Kit 1 Contains 11 items | 39 |
| S303 | Jobson Horne Probe | 40 |
| S304 | Ear Specula 4mm | 42 |
| S305 | Ear Specula 5mm | 42 |
| S306 | Ear Specula 6mm | 42 |
| S307 | Cerumen Hook | 40 |
| S308 | Rosen Inserter | 41 |
| S309 | Ear Specula 3mm | 42 |
| S310 | Myringotomy Knife | 39 |
| S311 | Jobson Horne Probe (Non-Sterile) | 40 |
| S312 | Cerumen Hook (Non-Sterile) | 40 |
| S313 | Ear Specula 3mm (Non Sterile) | 42 |
| S314 | Ear Specula 4mm (Non-Sterile) | 42 |
| S315 | Ear Specula 5mm (Non-Sterile) | 42 |
| S316 | Ear Specula 6mm (Non-Sterile) | 42 |
| S317 | Rosen Inserter (Non-Sterile) | 41 |
| S318 | Rounded Long Suction 12fg 20cm | 11 |
| S319 | Rounded Long Suction 9fg 20cm | 11 |
| S320 | Tischler-Morgan Biopsy Punch 200mm, Tip Down | 32 |
| S321 | Soft Tip 6fg 9cm Suction | 14 |
| S322 | Tumarkin Ear Specula Size 0.5, 4.5mm | 42 |
| S340 | Frazier Suction 15fg 9cm | 7 |
| S343 | Fenestrated Rounded Frazier 12fg 13cm | 10 |
| S353 | Pulse Lavage Kit includes Pulse Lavage and 2 tips | 43 |
| S367 | Kerrison Punch 2mm x 200mm @90° | 36 |
| S380 | Tilley Nasal Forceps Long | 38 |
| S419 | Yasargil Retractor | 34 |
| S432 | Insulated Suction 9fg 17cm 90° | 23 |
| S433 | Insulated Suction 6fg 17cm 90° | 23 |
| S434 | Micro Suction 18g 7cm | 16 |
| S435 | Retractor Suction 12fg 12cm 60° | 15 |
| S436 | Dissector Suction Flattened Tip 12fg 13cm | 15 |
| S438 | Diathermy Suction 9fg 17cm 70° | 25 |
| S439 | Teardrop Micro Suction 18g | 13 |
| S440 | Poole Tip Suction 12fg 13cm Rounded 4x4x2mm Fenestrations | 19 |
| S442 | Dissector Suction Flattened Tip 12fg 60° | 15 |
| S444 | Ball End Sinus Suction 9fg 60° Curve | 17 |
| S445 | Rounded Suction 12fg 30cm 4x1mm Fenestrations 60° | 22 |
| S446 | Curved Micro Suction 9fg 15cm 30° Bend 60° Curved Tip (Left) | 16 |
| S447 | Curved Micro Suction 9fg 15cm 30° Bend 60° Curved Tip (Right) | 16 |
| S448 | Diathermy Abbey Needle with Suction 6fg 10cm 40° Bend | 26 |
| S467 | Ball End Suction Upward Curve 70° 9fg | 17 |
| S468 | Ball End Suction Upward Curve 65° 12fg | 17 |




| Notes | Notes |
|-------|-------|
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |

| Notes | Notes |
|-------|-------|
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |



Other Product Ranges from Eakin Surgical

If you would like more detailed information on any of our other ranges, please contact your Eakin Surgical account manager or our customer services team who will be happy to help.

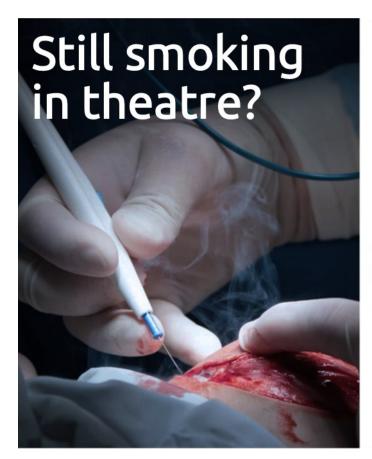



CIMPAX® C-VAC TORNADO™ Smoke Evacuation Electrosurgical Pencil

- Clears smoke as it is generated, improving air quality and visibility
- Range of accessories with ergonomic grip, ultrasoft push buttons and non-slip cover.











CONTACT US

Eakin Surgical Ltd Greypoint Cardiff Business Park Cardiff CF14 5WF United Kingdom

Telephone +44 (0)2920 767 800
Facsimile +44 (0)2920 767 801
Email info@eakinsurgical.com
Website www.eakinsurgical.com

Please send orders to orders@eakinsurgical.com

Please note products marked "*Made to Order" are not subject to standard delivery times and require a minimum order quantity. For more information, please contact info@eakinsurgical.com

These devices are intended for an individual patient, during a single procedure and should be discarded after use. Potential of cross contamination if re-used due to device not being sterile.

Eakin Surgical Ltd is registered in England & Wales. Company Registration No. 3939863. VAT No. GB806634432 Registered office: Eakin Surgical Ltd, Greypoint, Cardiff Business Park, Cardiff, CF14 5WF.

V14 10/2021

